

CUDC-101


#Cat: NB-64-19798-1ml Size: 1 ml #Cat: NB-64-19798-5mg Size: 5 mg #Cat: NB-64-19798-10mg Size: 10 mg #Cat: NB-64-19798-25mg Size: 25 mg #Cat: NB-64-19798-50mg Size: 50 mg #Cat: NB-64-19798-100mg Size: 100 mg Size: 200 mg #Cat: NB-64-19798-200mg #Cat: NB-64-19798-500mg Size: 500 mg

Chemical Properties:

CAS No: 1012054-59-9 **Formula:** $C_{24}H_{26}N_4O_4$ **Molecular Weight:** 434.49

Appearance: no data available

Storage: Powder: -20°C for 3 years | In solvent: -80°C for 1 year

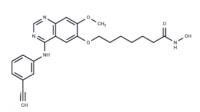

Biological Description:

Description	CUDC-101 is a potent inhibitor of HDAC, EGFR, and HER2 with IC50s of 4.4 nM, 2.4 nM,					
	and 15.7 nM, respectively					
Targets (IC50)	EGFR, HER, HDAC					
In vitro	CUDC-101 exhibits a dose-dependent inhibition of growth in erlotinib-sensitive H358 NSCLC xenografts and effectively suppresses tumor growth in erlotinib-resistant A549 NSCLC xenograft models. Additionally, CUDC-101 significantly promotes tumor regression in lapatinib-resistant, HER2-negative, EGFR-overexpressing MDA-MB-468 breast cancer models, and EGFR-overexpressing CAL-27 head and neck squamous cell carcinoma (HNSCC) models. Administration of CUDC-101 at a daily dose of 120 mg/kg in the Hep-G2 liver tumor model induces tumor deterioration more effectively than the maximum tolerated dose of erlotinib (25 mg/kg/day) and an equimolar dose of vorinostat (72 mg/kg/day).					
In vivo	CUDC-101 exhibits broad-spectrum antiproliferative activity in numerous human cancer cell lines with an IC50 range of 0.04 μ M to 0.80 μ M. It demonstrates higher potency compared to erlotinib, lapatinib, and combinations of these with either vorinostat or combinations of erlotinib or lapatinib. Notably, CUDC-101 can inhibit the proliferation of cancer cell lines resistant to lapatinib and erlotinib. This compound increases the acetylation of histones H3 and H4 and non-histone substrates like p53 and α -tubulin in a dose-dependent manner across various cancer cell lines. Additionally, CUDC-101 suppresses the expression of HER3, the amplification of Met, and the reactivation of AKT in tumor cells					



Kinase Assay	HDAC, EGFR and HER2 inhibition assays: The activities of Class I and II HDACs are assessed						
	using the Biomol Color de Lys system. Briefly, HeLa cell nuclear extracts are used as a source						
	of HDACs. Different concentrations of CUDC-101 are added to HeLa cell nuclear extracts						
	the presence of a colorimetric artificial substrate. Developer is added at the end of the						
	assay and enzyme activity is measured in the Wallac Victor II 1420 microplate reader at 405						
	nM. EGFR and HER2 kinase activity are measured using HTScan EGF receptor and Hi						
	kinase assay kits. Briefly, the GST-EGFR fusion protein is incubated with syntheti						
	biotinylated peptide substrate and varying concentrations of CUDC-101 in the presence of						
	400 mM ATP. Phosphorylated substrate is captured with strapavidin-coated 96-well plates.						
	The level of phosphorylation is monitored by antiphospho-tyrosine- and europium-labeled						
	secondary antibodies. The enhancement solution is added at the end of the assay a enzyme activity is measured in the Wallac Victor II 1420 microplate reader at 615 nM.						
Cell Research	Cancer cell lines are plated at 5000 to 10000 cells per well in 96-well flatbottomed plates						
	with varying concentrations of CUDC-101. The cells are incubated with CUDC-101 for 72						
	hours in the presence of 0.5% of fetal bovine serum. Growth inhibition is assessed by an						
	adenosine triphosphate (ATP) content assay using the Perkin-Elmer ATPlite kit. Apoptosis is						
	routinely assessed by measuring the activities of Caspase-3 and -7 using Apo-ONE						
	Homogeneous Assay Kit.(Only for Reference)						

Solubilit	У	Ethanol: 1 mg/mL (insoluble or slightly soluble),		
DMSO: 4.35 mg/mL (10 mM), Sonication is recommended.				
		(< 1 mg/ml refers to the product slightly soluble or insoluble)		

Preparing Stock Solutions

	1mg	5mg	10mg
1 mM	2.3015 mL	11.5077 mL	23.0155 mL
5 mM	0.4603 mL	2.3015 mL	4.6031 mL
10 mM	0.2302 mL	1.1508 mL	2.3015 mL
50 mM	0.046 mL	0.2302 mL	0.4603 mL

Please select the appropriate solvent to prepare the stock solution, according to the solubility of the product in different solvents. Please use it as soon as possible.

Reference

Cai X, et al. J Med Chem, 2010, 53(5), 2000-2009. Lai CJ, et al. Cancer Res, 2010, 70(9), 3647-3656.

> $Inhibitor \cdot \textbf{Natural Compounds} \cdot \textbf{Compound Libraries} \cdot \textbf{Recombinant Proteins}$ This product is for Research Use Only · Not for Human or Veterinary or Therapeutic Use